

OPEN SURGERY



SWING SYSTEM Laminectomy Punches

WORLDWIDE patented One-piece detachable Punches

HIGH PUNCHING PERFORMANCE

Extraordinary precision due to perfect and accurate workmanship

Easy removal of vertebral pieces during surgery

BEFORE & AFTER SURGERY: EASY HANDLING

Quick and intuitive assembly and disassembly

The upper slider is hinged: No mix-up or loose of slider parts during cleaning and sterilization

EFFICIENT CLEANING

The take-apart design allows to clean the inner and outer surfaces of the instrument

Thorough cleaning, disinfection and sterilization prolongs the service life of the instrument.





SWING SYSTEM

Laminectomy Punches

WORLDWIDE patented One-piece detachable Punches

All SWING Laminectomy punches are available as Standard, Thin Footplate ("T") and with Ceramic Coating ("CC"). Punches with 2-6 mm width of bite can be ordered with or without ejector ("A"), only the punches with 1 mm width of bite are available only without ejector.

Example how to order punches with eligible options (ejector, Thin footplate, ceramic coated):

You would like to order the 249-071-02 with ejector ("A"), thin footplate ("T") and ceramic coated ("CC")?


► Order: number 249-071-02ATCC

INTUITIVE HANDLING


Assembly



First push shaft parts into one other ...

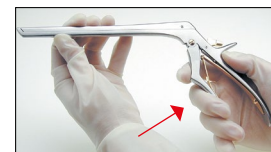

then close the handle ...

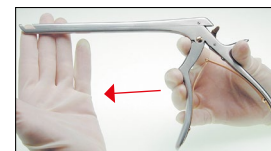


and finally open the handle

Disassembly



First close the handle (1); then push the lever downwards (2)...

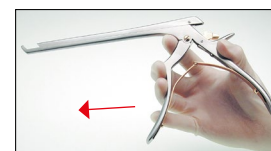

now open the handle ...



and push downwards at the end of the shaft

100 % Made in Germany

KERRISON SWING SYSTEM - Upcut 90°; 10 mm open

The ejector is included for instruments with width of bite 2-6 mm

Width of bite	Shaft length			Ejector	Thin footplate	Ceramic coated
	18.0 cm; 7.2"	20.0 cm; 8.0"	23.0 cm; 9.2"			
1 mm	249-051-01	249-052-01	249-053-01	none	T	CC
2 mm	249-051-02	249-052-02	249-053-02	incl.	T	CC
3 mm	249-051-03	249-052-03	249-053-03	incl.	T	CC
4 mm	249-051-04	249-052-04	249-053-04	incl.	T	CC
5 mm	249-051-05	249-052-05	249-053-05	incl.	T	CC
6 mm	249-051-06	249-052-06	249-053-06	incl.	T	CC

KERRISON SWING SYSTEM - Upcut 40°; 10 mm open

The ejector can be selected for instruments with width of bite 2-6 mm

Width of bite	Shaft length			Ejector	Thin footplate	Ceramic coated
	18.0 cm; 7.2"	20.0 cm; 8.0"	23.0 cm; 9.2"			
1 mm	249-071-01	249-072-01	249-073-01	none	T	CC
2 mm	249-071-02	249-072-02	249-073-02	A	T	CC
3 mm	249-071-03	249-072-03	249-073-03	A	T	CC
4 mm	249-071-04	249-072-04	249-073-04	A	T	CC
5 mm	249-071-05	249-072-05	249-073-05	A	T	CC
6 mm	249-071-06	249-072-06	249-073-06	A	T	CC

KERRISON SWING SYSTEM MODIFIED - Upcut 40°; 15 mm open

Without ejector; thin footplate; ceramic coating optional

Width of bite	Shaft length	
	20.0 cm; 8.0"	add "CC"
2 mm	249-272-02T	CC
3 mm	249-272-03T	CC
4 mm	249-272-04T	CC
5 mm	249-272-05T	CC



TONTARRA Medizintechnik GmbH
 Daimlerstr. 15 · 78573 Wurmlingen/Germany
 Tel: +49(0) 74 61/9 65 76 -0) · info@tontarra.de
www.tontarra.de

Authorized distributor for

Distributed in the UK by



7-9 Parkside Centre, Potters Way,
 Southend-on-Sea
 Essex SS2 5SJ T 01702
 602050 E
info@surgicalholdings.co.uk

All rights are reserved. Technical alterations are possible. No part may be copied or reproduced in any form. In the case of misuse, we retain the rights to recall our promotional materials and products and to take legal actions. 999-115-02E 06.2018